

**DOUBLE
BLACK
IMAGING**


5MP Mono LCD

DBI2105MP-PP

The DBI2105MP-PP is a 5 mega-pixel high bright monochrome LCD offering fixed front sensors for True DICOM gray level and luminance auto-calibration. The unobtrusive sensor emerges from the bezel to calibrate and retracts when calibration is complete.

X-CAL® Client and Remote Administration DICOM Auto-Calibration Software are included with each system at no charge. X-CAL® Client, the local workstation software will allow even individual workstations to be configured for auto-calibration without a network software piece. X-CAL® Remote Administration Server Software and the P-Series will allow an administrator to monitor all the displays within the enterprise, launch calibration and receive alerts if LCDs fall out of calibration or if calibration was unable to be completed.

The durable, bonded protective panel includes an anti-reflective coating and comes standard with the DBI2105MP-PP. This added durability reduces maintenance costs and extends the useful life of the LCD.

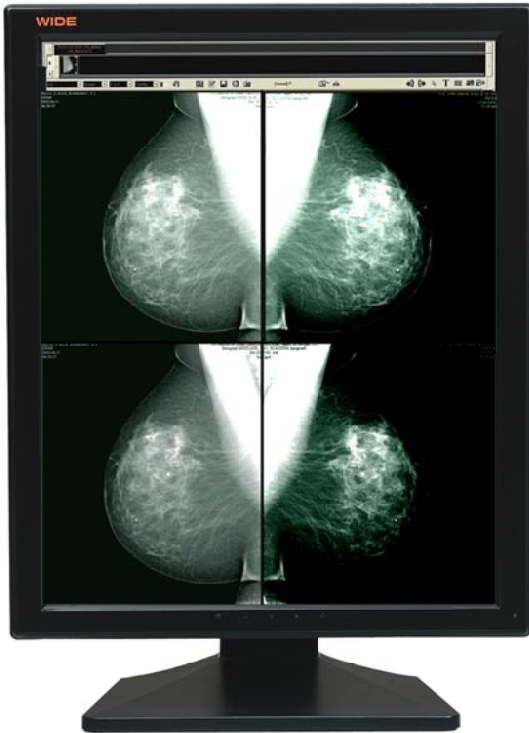

- ◆ 21.3" 5MP (2560 x 2048, 2048 x 2560) Grayscale TFT LCD
- ◆ Brightness - 1000 cd/m² (Typical); Contrast Ratio - 600:1 (Typical)
- ◆ Complete Screen Luminance Uniformity Control
- ◆ Built in retractable front sensors with backlight sensor
- ◆ 12 Bit Processing for Precise Image Rendering
- ◆ Digital Ambient Light Compensation
- ◆ Bonded Protective Panel
- ◆ Full 5 Year Warranty

Model No. DBI2105MP-PP

Description	5 Mega-Pixel LCD with Fixed Front Sensors, Backlight Sensor, and Bonded Protective Panel
Screen Type	Flat Panel LCD
Recommended Modalities	Mammography, CR, CT, DR, MR, US, PACS
Diagonal Screen Size (inches)	21.3"
# of Pixels	5.2 Mega-Pixel
Resolution (pixels)	2048 x 2560 Portrait; 2560 x 2048 Landscape
Contrast Ratio (typical)	600:1
Viewing Angle	170° horizontal & vertical
Brightness	1000 cd/m ²
Palette	12 Bit
Video Inputs	DVI
Dimensions H x W x D (inches)	21.6 x 15.6 x 7.8
Weight with stand/base (lbs.)	22.7
Power Supply	External for added reliability
Power Consumption	115W
Compliance/Safety Standards	UL60601-1, cUL2601-1, CE, IEC/EN60601-1, FCC Class B, VCC1, MIC, CCC, FDA 510(k)
Service/Support	24/7
Warranty	5 Years on the LCD, Electronics and Back light; 3 years on the graphic controller
Other Features	Retractable Front Sensor for True Auto-Calibration, Client/Remote Administration Calibration Software included, Bonded Protective Panel, Rear Sensor for Back light brightness control, Digital Ambient Light Compensation, Complete Screen Luminance Uniformity Control, 12 Bit Processing for Precise Image Rendering, VESA Standard Mounting configuration, Other Mounting options available



Double Black Imaging

11005 Dover St., Suite 500
Westminster, CO 80021
Tel. 877-852-2870, 303-404-2222
Fax 303-460-0276
sales@doubleblackimaging.com
www.doubleblackimaging.com

