


X-RAY GENERATORS

CMP 200[®] DR

(SERIES)

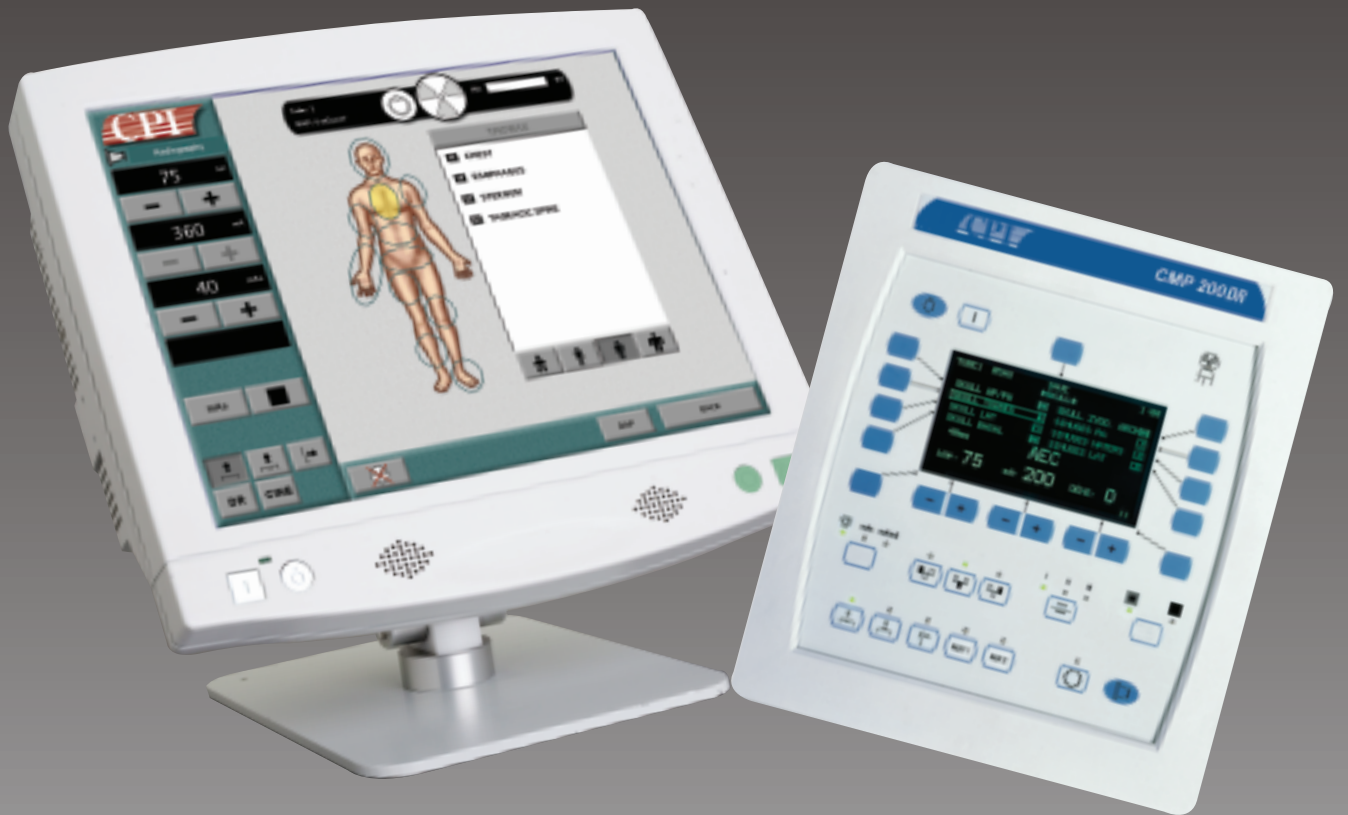

CMP 200[®] DR series

Radiographic 200 kHz⁺ X-Ray generators – unrivalled reliability and performance.

4 different models to choose from:

- 40 kW
 - 50 kW
 - 65 kW
 - 80 kW
- Digital Radiography Interfaces*
 - Dual Speed Starter*
 - Cost-effective for all radiographic applications (hospitals, medical clinics, chiropractic and veterinary applications)
 - Field upgradable for your future imaging requirements.

**Small. Compact. Lightweight.
Cost Effective.**

* Consult factory for details and other options

X-RAY
**Generators
of Choice**

The newest n CPI family 'Generators



Membrane Console
Standard

'Consoles-

member to the of X-Ray s-of-Choice'.



Touch Screen Optional

of-Choice'

Look at these Features

- CMP200® DR generator utilizes CPI's proven 200 kHz+ inverter technology and software architecture.
- Anatomical programs available in many languages that are easily edited by the Operator Control Console
- Technique Selection:
 1. kV with AEC*
 2. kV/mAs
 3. kV/mA/Time
 4. Patient/Body Part Thickness (Membrane Control Console)
- Low kV ripple – maximizes image quality by providing constant potential performance
- Automatic Tube Calibration
- GenWare® Service Software allows fast and easy set up as well as system diagnostics for off-site analysis
- Collimator Power Supply
- Basic Tomography Interface

Options:

- Digital Radiography Interfaces*
- Touch Screen Control Console (with embedded GenWare® Service Software)
- Dual Speed Starter for higher power
- Exposure Hand Switch
- Mini Console (for Integrated Systems)
- Console Wall Mount
- Automatic Exposure Control (AEC) Interface
- Dose Area Product (DAP) Interface
- Auto Positioning Interface

Compliance with FDA, UL/CSA, MDD, CE, CQC, and SFDA Requirements

* Consult factory for details and other options



Exposure Hand Switch Optional



Mini Console Optional



Wall Mount
Optional



Membrane Console Pedestal
Optional

World Class Product Support

One of our greatest strengths is our world class reputation for customer support. A dedicated support team is only a phone call away at 905-877-0161 – or at our North American toll-free number 1-888-CPI-XRAY (274-9729)

Call us. We are here to help you with all of your installation needs.



Communications & Power Industries Canada Inc.
communications & medical products division
45 River Drive / Georgetown, Ontario / Canada L7G 2J4
Hot Line Telephone: 1-888-CPI-XRAY
TEL: 905-877-0161 / FAX: 905-877-5327
E-MAIL: marketing@cmp.cpii.com
WEB: www.cpii.com



Some features shown may be optional. Please refer to Data Sheet.